


# Validation Data Sheet

QuantiNova LNA PCR Assay	KRT14
Assay Catalog No.	SBH1220158
dPCR Instrument	QIAcuity digital PCR instrument
dPCR Master Mix	QIAcuity EG PCR Master Mix
Experimental Sample(s)	gBlock

**Gene Expression Analysis:** Concentration values [copies/ $\mu$ l] for different gBlock sample amounts (0 copies, 4400 copies, 6600 copies, 8880 copies and 13320 copies). Technical triplicates are shown.

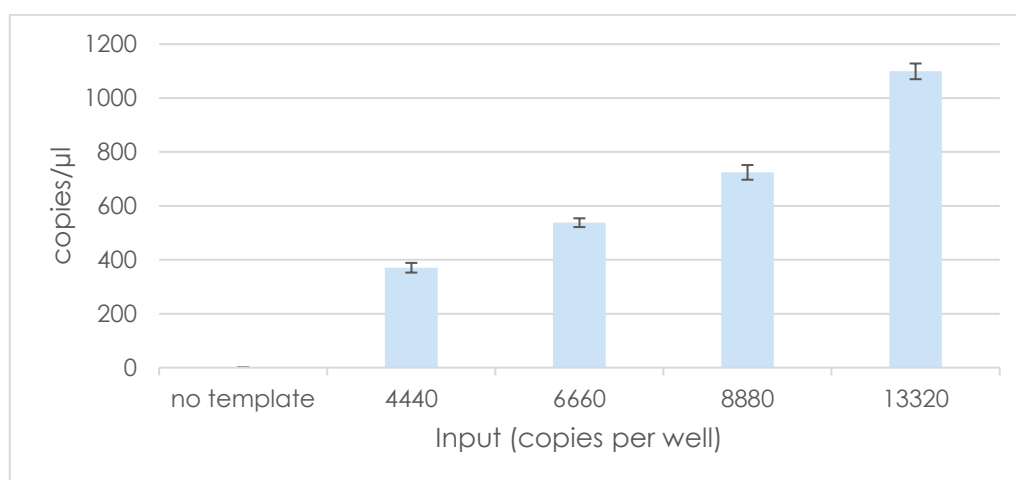

**1D Scatter Plot:** Single-well data (representative well of technical triplicates) for the different gBlock input amounts (FAM channel for EvaGreen detection; gBlock amounts in the 5 wells: 0 copies/ $\mu$ l, 370 copies/ $\mu$ l, 555 copies/ $\mu$ l, 740 copies/ $\mu$ l and 1110 copies/ $\mu$ l).

